



Yes, you can.

SC900DLX



SHOWN WITH OPTIONAL PAN DECK


The SC900DLX is a full-function, full-feature, low-height bed that promotes quality resident care with full bed performance. Its three DC motors provide ultra-quiet operation of the head, hi-low and foot functions of the bed. Plus, the horizontal travel range of the SC900DLX is designed so that the bed will not damage walls during height adjustments.

SC900DLX



6-Function Hand Control Pendant

The electric Hand Control Pendant features bright, universal symbols and allows residents and staff to easily operate the comfort controls of the bed.

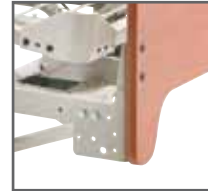

Adjustable Mattress Deck

With adjustable bed-end brackets and corner mattress keepers, the SC900DLX beds are expandable from 76" to 80".

SC900DLX Specifications

- Deck length: 76" or 80"
- Overall length: 82" with bed ends
- Sleep surface deck to floor- low position: 9⁷/₈"
- Sleep surface deck to floor- high position: 20¹/₄"
- Under-bed clearance- high position: 5"
- Overall width (with rails): 36"
- Weight capacity: 450 lbs
- Emergency Override: In the event of a power outage, or if the battery runs out, the SC900DLX features a built-in hand crank so it can be raised or lowered by hand.

Accessories

- Vinyl bumpers
- Trapeze/mounting frame
- Pan Deck

Invacare Continuing Care, Inc.

www.invacare-ccg.com

info_iccg@invacare.com

USA

One Invacare Way
Elyria, Ohio 44035
USA
800-668-2337

Canada

994 Hargrievie Road
London, Ontario
N6E 1P5 Canada
800-668-2337

Trademarks are identified by the symbol™ and ®.
All trademarks are owned or licensed to Invacare Corporation or its subsidiaries unless otherwise noted.

© 2009 All rights reserved. Invacare Corporation.

Form 08.40.9/09



Yes, you can: